



Refrigerator Instructions for Use

i.Series® · Horizon Series™
Upright - Undercounter - Pass-Thru

Blood Bank

i.Series

iBR105-GX, iBR113-GX, iBR120-GX,
iBR125-GX, iBR245-GX, iBR256-GX,
iBR226-GX, iBR458-GX

Horizon Series

HBR105-GX, HBR113-GX, HBR120-GX,
HBR125-GX, HBR245-GX, HBR256-GX,
HBR226-GX, HBR458-GX

Laboratory

i.Series

iLR105-GX, iLR113-GX, iLR120-GX,
iLR125-GX, iLR245-GX, iLR256-GX

Horizon Series

HLR105-GX, HLR113-GX, HLR120-GX,
HLR125-GX, HLR245-GX, HLR256-GX

Pharmacy

i.Series

iPR105-GX, iPR113-GX, iPR120-GX,
iPR125-GX, iPR245-GX, iPR256-GX,
iPR226-GX, iPR458-GX

Horizon Series

HPR105-GX, HPR113-GX, HPR120-GX,
HPR125-GX, HPR245-GX, HPR256-GX,
HPR226-GX, HPR458-GX



Document History

Revision	Date	CO	Supersession	Revision Description
A	18 SEP 2019*	14979	n/a	Initial release.
B	19 MAY 2020*	15365	B supersedes A	<ul style="list-style-type: none">Updated Location Requirements to include reference to Product Specifications.Updated Chart Recorder section to include reference to Temperature Chart Recorder manual.Updated Notified Body information in the Compliance section.
C	25 JAN 2022*	16229	C supersedes B	<ul style="list-style-type: none">Updated exterior dimensions in Refrigerator Dimensions table.
D	12 OCT 2022*	17064	D supersedes C	<ul style="list-style-type: none">Updated content to include Pass-Thru Refrigerator information
E	9 DEC 2025	30216	E supersedes D	<ul style="list-style-type: none">Updated the Compliance section to reflect the most current requirements.

* Date submitted for Change Order review. Actual release date may vary.

Document Updates

The document is furnished for information use only, is subject to change without notice and should not be construed as a commitment by Helmer Scientific. Helmer Scientific assumes no responsibility or liability for any errors or inaccuracies that may appear in the informational content contained in this material. For the purpose of clarity, Helmer Scientific considers only the most recent revision of this document to be valid.

Notices and Disclaimers

Confidential / Proprietary Notices

Use of any portion(s) of this document to copy, translate, disassemble or decompile, or create or attempt to create by reverse engineering or otherwise replicate the information from Helmer Scientific products is expressly prohibited.

Copyright and Trademark

Copyright © 2025 Helmer, Inc. Helmer®, i.Series®, i.C³®, Horizon Series™, and Rel.i™ are registered trademarks or trademarks of Helmer, Inc. in the United States of America. All other trademarks and registered trademarks are the property of their respective owners. Helmer, Inc., doing business as (DBA) Helmer Scientific and Helmer.

Disclaimer

This manual is intended as a guide to provide the operator with necessary instructions on the proper use and maintenance of certain Helmer Scientific products.

Any failure to follow the instructions as described could result in impaired product function, injury to the operator or others, or void applicable product warranties. Helmer Scientific accepts no responsibility or liability resulting from improper use or maintenance of its products.

The screenshots and component images appearing in this guide are provided for illustrative purposes only, and may vary slightly from the actual software screens and/or product components.

Contents

- 1 About this Manual 1**
 - 1.1 Intended Audience 1
 - 1.2 Model Reference 1
 - 1.3 Intended Use 1
 - 1.4 Safety Precautions and Symbols 1
 - 1.5 Avoiding Injury 3
 - 1.6 General Recommendations 3
- 2. Installation 4**
 - 2.1 Location 4
 - 2.2 Placement and Leveling 4
 - 2.3 Stacked Undercounter Units 4
 - 2.4 AC Power Cord 5
 - 2.5 Temperature Probes 5
 - 2.6 Chart Recorder (if included) 6
- 3 i.Series® Operation 8**
 - 3.1 Initial Power-Up 8
 - 3.2 Operation 9
 - 3.3 Change Temperature Setpoint 9
 - 3.4 Set Alarm Parameters 10
 - 3.5 Active Alarms 11
 - 3.6 Mute Active Alarms 11
 - 3.7 Light Operation (if installed) 12
- 4 Min/Max Temperature Monitoring 13**
- 5 i.Series® Access Control and i.Lock™ Interlock Feature (options for select models) 14**
 - 5.1 Using Access Control (not available on Pass-Thru) 14
 - 5.2 Using i.Lock™ Interlock (Pharmacy Pass-Thru only) 15
- 6 Horizon Series™ Operation 16**
 - 6.1 Initial Power-Up 16
 - 6.2 Display Minimum and Maximum Monitor Temperature Recordings 16
 - 6.3 Change Temperature Setpoint 17
 - 6.4 Set Parameter Values 18
 - 6.5 Set Temperature Units 18
 - 6.6 Set Alarm Setpoints (Parameters) 18
 - 6.7 Temperature Calibration Offsets 19
 - 6.8 Active Alarms 19
 - 6.9 Mute and Disable Audible Alarms 19
 - 6.10 Light Operation (if installed) 19
- 7 Horizon Series™ Access Control (Optional) 20**
 - 7.1 Using Access Control 20
- 8 Product Specifications 21**
 - 8.1 Operating Standards 21
 - 8.2 Refrigerator Specifications 22
- 9 Compliance 22**
 - 9.1 Safety Compliance 22
 - 9.2 Environmental Compliance 22
 - 9.3 EMC Compliance 22

1 About this Manual

1.1 Intended Audience

This manual provides information on how to use i.Series® and Horizon Series™ refrigerators. It is intended for use by end users of the refrigerator and authorized service technicians.

1.2 Model Reference

Models are indicated by a distinguishing model number that corresponds to the series, type, number of doors, and capacity of the refrigerator. For example, “iLR113-GX” refers to an i.Series Laboratory Refrigerator with one door and a capacity of 13 cu ft and “HPR458” refers to a Horizon Series Pharmacy Pass-Thru Refrigerator with four doors and a capacity of 58 cu ft.

1.3 Intended Use

Note

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at their own expense.

Helmer refrigerators are intended for the storage of blood products and other medical and scientific products.

1.4 Safety Precautions and Symbols

Symbols found in this document

The following symbols are used in this manual to emphasize certain details for the user:



Task Indicates procedures which need to be followed.



Note Provides useful information regarding a procedure or operating technique when using Helmer Scientific products.

NOTICE Advises the user against initiating an action or creating a situation which could result in damage to equipment; personal injury is unlikely.

Symbols found on the units

The following symbols may be found on the refrigerator or refrigerator packaging:



Caution: Risk of damage to equipment or danger to operator



Danger: Risk of Fire or Explosion. Flammable refrigerant used



Caution: Hot surface



Refer to documentation



Caution: Shock / electrical hazard



EU Authorized Representative



Warning: Crushing of hands / fingers

These symbols also appear with appropriate information provided within this document.

1.5 Avoiding Injury



- Do not use mechanical devices or other means to accelerate the defrosting process, other than those recommended by the manufacturer.
- Do not damage the refrigerant circuit.

Review safety instructions before installing, using, or maintaining the equipment.

- ◆ Before moving unit, ensure door is closed and casters are unlocked and free of debris.
- ◆ Before moving unit, disconnect the AC power cord and secure the cord.
- ◆ Pass-Thru units - Do not move a unit if load exceeds 900 lbs/408 kg (226 units) or 1350 lbs/612 kg (458 units).
- ◆ Never physically restrict any moving component.
- ◆ Avoid removing electrical service panels and access panels unless so instructed.
- ◆ Keep hands away from pinch points when closing doors.
- ◆ Avoid sharp edges when working inside the electrical compartment and refrigeration compartment.
- ◆ Ensure products are stored at recommended temperatures determined by standards, literature, or good laboratory practices.
- ◆ Proceed with caution when adding and removing product from the refrigerator.
- ◆ Do not open multiple, loaded drawers at the same time.
- ◆ Use manufacturer supplied power cord only.
- ◆ Avoid risk of ignition by using only manufacturer supplied components and authorized personnel when servicing the unit.
- ◆ Using the equipment in a manner not specified by the manufacturer may impair the protection provided by the equipment.
- ◆ Ensure product is stored safely, in accordance with all applicable organizational, regulatory and legal requirements.
- ◆ The refrigerator is not considered to be a storage cabinet for flammable or hazardous materials.
- ◆ Use caution when moving undercounter units in a stacked configuration.
- ◆ **REQUIRED:** Decontaminate parts prior to sending for service or repair. Contact Helmer or your distributor for decontamination instructions and a Return Authorization Number.

1.6 General Recommendations

General Use

Allow refrigerator to come to room temperature before switching power on.

During initial startup, high temperature alarm may sound while refrigerator reaches operating temperature.



Do not remove the cover from the condensate evaporator tray on upright units.

Initial Loading

Allow chamber temperature to stabilize at the setpoint before storing product.

Product Loading Guidelines

When loading your refrigerator, take care to observe the following guidelines:

- ◆ Never load refrigerators beyond capacity.
- ◆ Always store items within shelves or drawers.
- ◆ Temperature uniformity is maintained by air circulation, which could be impeded if unit is overfilled, particularly at the top or against the doors or walls. Ensure a 2" (50 mm) clearance is provided below the fan.

Note

Products stacked against walls or doors may obstruct air flow and affect performance of unit.

2. Installation

2.1 Location



Keep all ventilation openings in the enclosure or, in the structure of building-in, clear of obstruction.

- ◆ Has a grounded outlet meeting the electrical requirements listed on the product specification label.
- ◆ Is clear of direct sunlight, high temperature sources, and heating and air conditioning vents.
- ◆ Meets specified limits for ambient temperature and relative humidity as stated in the Product Specifications section of this manual.
- ◆ Upright units require minimum 8" (203 mm) above, and minimum 3" (76 mm) behind.
- ◆ Undercounter units require minimum 3" behind the unit for clearance and feature access.
- ◆ Pass-Thru units - Clearance above the clean room side may be 0"; minimum clearance above the ante-room side is 8" (203 mm) .
- ◆ Pass-Thru units - Minimum clearance between the cabinet and wall should be 0".
- ◆ Pass-Thru units should be positioned so the wall header is over the solid cover, but not over the flange or blue bezel.
- ◆ Pass-Thru units may need additional clearance for the following scenarios: placing a unit where hinges are against the wall; placing units side-by-side; using an optional trim kit on a single unit, or on one or both adjacent units.
- ◆ Pass-Thru units - Clearance needed between adjacent units:
 - 6.2" for full access to both units with trim around each unit
 - 5.0" for full access to both units with trim around one unit
 - 2.5" for full access to one unit with no trim installed.

2.2 Placement and Leveling



The evaporation tray located on the back of the upright refrigerator may be hot. Do not use the tray as a handle.

NOTICE

- Helmer does not recommend operating this unit on a GFI/GFCI outlet.
- To prevent tipping, ensure the casters (if installed) are unlocked and doors are closed before moving the unit.
- To avoid damaging refrigerant tubing or risking refrigerant leak, use caution when moving or operating the unit.

Undercounter units only

- Do not sit, lean, push or place heavy objects on top surface of undercounter units.
- Do not lean on or push down on an open door or extended drawers.

1. Ensure door(s) is secured and casters (if installed) are unlocked.
2. Roll refrigerator into place and lock casters.
3. Ensure refrigerator is level.

Note

Helmer recommends the use of leveling feet and wall and floor brackets for stabilization on undercounter units. Contact Helmer Technical Service for parts and instruction.

2.3 Stacked Undercounter Units

NOTICE

- For stacked configuration, both units must have leveling feet installed.
- Stacking kit includes back brace bars and front stabilizing brackets. All pieces must be installed.
- When stacking units, place the heavier unit on the bottom.
- Do not open multiple loaded drawers at the same time.
- Do not lean on or push down on an open door or extended drawers.

Contact Helmer or your distributor for more information regarding the stacking kit and methods to secure both units to the wall and/or floor.

2.4 AC Power Cord



Use manufacturer supplied power cord only.

Install power cord

If packaged with modular cord, insert plug securely into the refrigerator power receptacle prior to connecting to grounded outlet.

2.5 Temperature Probes

A solid ballast or probe bottle and container of glycerin have been provided with this unit. The glycerin is used to create a solution which, when placed in the probe bottle, simulates the product stored in the refrigerator. The product simulation solution temperature reflects the product's temperature during normal operation.

Notes

- Temperature probes are fragile; handle with care.
- Number and location of probes varies by model.
- Remote probes may also be introduced through the existing port on top or rear of the unit (if included).
- Solid ballast (if installed) should be placed in the bracket in a horizontal position.
- Failure to fill probe bottles or keep probe bottles (if installed) filled to the appropriate level may cause the chamber temperature to display higher or lower than the actual temperature.

Primary Monitor Probe

The primary monitor probe is located at the top left side of the refrigerator in single-door upright and undercounter units, and in the center, between doors, of the double-door units.

The primary monitor probe in pass-thru units, as accessed from the ante-room side, is located at the top left side of the single-door (226) units and in the center, between doors, of the double-door (458) units.



Primary monitor probe locations

Secondary Monitor Probe (*i.Series models 20 cu ft and larger only*)

The secondary monitor probe is located in the lower left side of the refrigerator in single-door and double-door upright units.

The secondary monitor probe in pass-thru units, as accessed from the ante-room side, is located in the lower right side of the single-door (226) and double-door (458) units.



Secondary monitor probe

☑ Fill Temperature Probe Bottle (if installed)**i Note**

Use approximately 4 oz. (120 mL) of product simulation solution (10:1 ratio of water to glycerin). Glycerin packet included in refrigerator box.

1. Remove probe(s) from bottle and remove bottle from bracket.
2. Remove cap and fill with approximately 4 oz. (120 mL) of product simulation solution.
3. Secure cap on bottle and place in bracket.
4. Replace probe(s), immersing at least 2" (50 mm) in solution.

☑ Install Additional Probe Through Access Port

1. Peel back putty to expose port.
2. Insert probe through port into chamber.
3. Insert probe into bottle.
4. Replace putty, ensuring a tight seal.

2.6 Chart Recorder (if included)

The chart recorder has a back-up battery system enabling a period of continuous operation if power is lost. Battery life varies by manufacturer as well as voltage level remaining. If full battery power is available, back-up power for the temperature chart recorder is available for up to 14 hours.

i Notes

- If chart recorder is operated on battery power, the battery should be replaced to ensure the back-up source has proper charge.
- For complete information, refer to the Temperature Chart Recorder Operation and Service Manual.

Prior to use:

Place chart recorder probe in bottle or ballast with primary monitor probe.

Set up and Operation

Access the chart recorder by pressing and releasing the door (i.Series except 113 models) or pulling door open (Horizon Series, 113 models and undercounter models).



☑ Install battery.

Connect the leads to the battery to provide back-up power to the chart recorder.

☑ Install / Replace Chart Paper**i Note**

For accurate temperature reading, ensure the current time is aligned with the time line groove when the chart knob is fully tightened.



Chart recorder stylus and time line groove

1. Press and hold C button. When stylus begins to move left, release button. The LED flashes.
2. When stylus stops moving, remove chart knob then move knob up and away.
3. Place chart paper on chart recorder.
4. Gently lift stylus and rotate paper so current time line corresponds to time line groove.
5. Hold chart paper in place while making sure the chart knob is fully tightened. *(Failure to fully tighten the knob can result in paper slipping and losing time.)*
6. Press and hold C button. When stylus begins to move right, release button.
7. Confirm stylus is marking on paper and stops at the correct temperature.
8. Calibrate chart recorder to match primary temperature if needed and close recorder door.

Power Supply

The temperature chart recorder uses AC power when the system is operating. If AC power fails, the recorder continues to record temperature with back-up power provided by the nine-volt battery.

- ◆ The LED indicator glows green continually when main power is functioning and the battery is charged.
- ◆ The LED indicator glows red continually when main power is functioning and the battery is either not installed or needs to be changed.
- ◆ The LED indicator flashes red to indicate that the recorder is receiving power only from the back-up battery.
- ◆ The LED indicator flashes during chart paper change mode.

3 i.Series® Operation

3.1 Initial Power-Up

1. Plug the power cord into a grounded outlet that meets the electrical requirements on the product specification label.
2. Turn the AC power switch ON.
3. Turn the back-up battery switch ON.


i Notes

- For models equipped with optional Access Control, the back-up battery is turned ON with a key switch.
- The Start screen is displayed when the i.C³ is powered on. The i.C³ will take approximately 2-5 minutes to boot up.



Start screen

The language screen is displayed when the i.C³ is powered on. Use the Language screen to select the i.C³ display language.

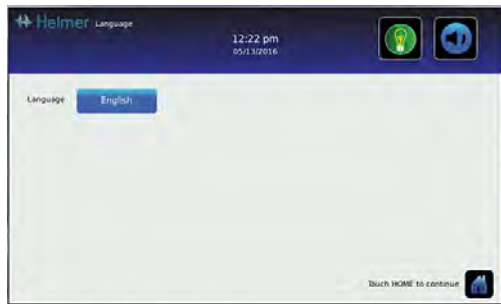

Language screen

If an alarm sounds, temporarily mute the alarm by touching the Mute button.



Home screen - alarm muted



Mute icon

i Note

Active alarms are displayed on the Home screen. If an alarm condition other than High Temperature occurs, refer to the service manual for troubleshooting.




3.2 Operation

Notes

- Refer to the i.C³ User Guide for complete information regarding the i.C³ User Interface.
- The i.C³ Home screen displays temperature and alarm information, and provides icons to gain access to other functions of the i.C³.
- After two minutes of inactivity, the screensaver will be displayed. To return to the Home screen, touch the screensaver.



Home screen

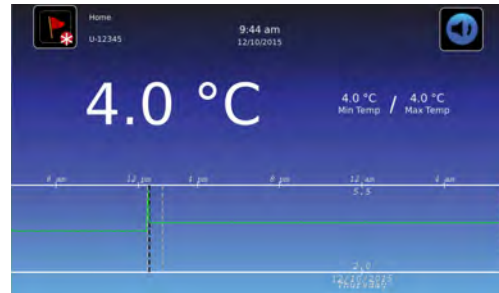

Home screensaver


3.3 Change Temperature Setpoint

Note


The Temperature Setpoint toggle button can be accessed from either the initial Settings screen or the Device Control Settings screen.



> Enter the Settings password. Select Temperature Setpoints. Touch minus (-) or plus (+) on spin box to change value.

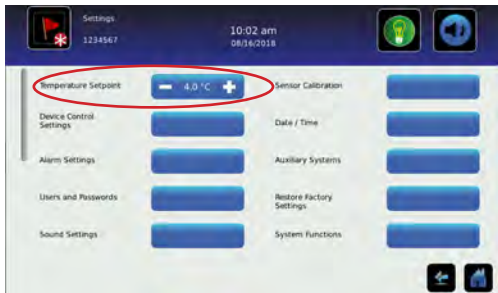

Settings screen

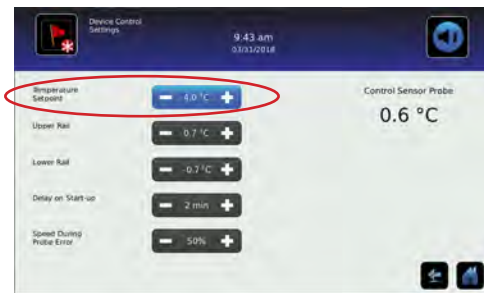



Device Control Settings screen

Notes

- Default Settings password is 1234.
- Default setpoint is 4.0 °C for iLR and iBR models, or 5.0 °C for iPR models.

3.4 Set Alarm Parameters

 >  > Enter the Settings password. Scroll down to select Alarm Settings. Touch minus (-) or plus (+) on the spin box to set each alarm parameter.



Settings screen



Alarms screens

Alarm settings control the circumstances and timing of alarm condition indicators displayed on the i.C³ Home screen.

i.Series

3.5 Active Alarms



Home screen with active alarm

Table 1. i.Series Active Alarms

Alarm	Description
Communication Failure 1	Communication lost between i.C ³ display board and control board
Communication Failure 2	Configuration file is corrupt or i.C ³ is unable to access the configuration file
Communication Failure 3	Corrupt database
Compressor Probe Failure	Probe not functioning properly
Compressor High Temperature	Compressor temperature reading is above high temperature alarm setpoint
Control Probe Failure	Probe not functioning properly
Drive Space Low	SD card is approaching capacity
Drive Space Full	SD card is full
Door Open	Door is open beyond user-specified duration
Inverter Communication Failure	Communication is lost between the i.C ³ control board and the VCC inverter
Low Battery	Back-up battery voltage is low
No Battery	Back-up battery voltage is deficient
Power Failure	Power to the unit has been disrupted
Primary Monitor Probe Failure	Probe not functioning properly
Primary Probe High Temperature	Primary monitor probe reading is above high temperature alarm setpoint
Primary Probe Low Temperature	Primary monitor probe reading is below low temperature alarm setpoint
Secondary Monitor Probe Failure (if installed)	Probe not functioning properly
Secondary Probe High Temperature (if installed)	Secondary monitor probe reading is above high temperature alarm setpoint
Secondary Probe Low Temperature (if installed)	Secondary monitor probe reading is below low temperature alarm setpoint

3.6 Mute Active Alarms

Audible alarms may be temporarily muted by touching the Mute icon. The delay duration can be set and changed by selecting Sound Settings from the Settings screen. The duration may be set to any value from 1 - 60 minutes. The delay time remaining will be displayed in the bottom right corner of the icon. If the alarm is still active after the mute delay has ended, the audible alarm will sound.



Unmuted

Muted



> Enter the Settings password. Scroll down to select Sound Settings. Touch minus (-) or plus (+) on spin box to set the mute duration.

3.7 Light Operation (if installed)

Press Light Icon to turn LED lights ON or OFF.
Auto ON/OFF feature can be configured in Settings.



Light ON/OFF

Table 2. Application Icons

Icon	Description	Icon	Description	Icon	Description	Icon	Description
	Home		Temperature Graph		CSV Download		Save
	Event Log		Alarm Test		PDF Download		Cancel
	Mute		Information Logs		Upload		Back Arrow
	Reset		Contact Information/ Contact Helmer		Access Control		Scroll
	Zoom Information		Display Brightness		Access Log		Temperature Graph Forward/Back
	i.C³ Applications		Light On/Off		Alarm Conditions		Battery Power
	Settings		Icon Transfer		Cancel Test		



4 Min/Max Temperature Monitoring

The Min/Max temperature display provides the highest and lowest Primary Monitor probe temperature reading since the last system reset (power-on event) or manually-initiated reset. Touch the Reset icon to the right of the display to manually reset.



Notes

- The Min/Max temperature display can be turned on or off through Display Settings.
- Once the time reaches the maximum display of 999 hours and 60 minutes, the message will display ">999:60", but minimum and maximum temperatures will continue to be tracked.

5 i.Series® Access Control and i.Lock™ Interlock Feature (options for select models)

Access Control and Interlock features limit access to the refrigerator.

5.1 Using Access Control (not available on Pass-Thru)

Notes

- During a power failure, the optional Access Control lock will remain locked until battery power is depleted or until the back-up battery key switch is switched OFF.
- Switching the back-up battery key switch OFF will disable the monitoring system during a power failure.
- During a power failure, switch the battery back-up switch OFF and use the mechanical door key to provide secure storage for refrigerator contents.
- Refer to i.C³ User Guide for complete information regarding Access Control.

Configure and manage user-specific accounts to allow controlled access to the refrigerator. Enter the supervisor PIN to set up Access Control and follow the on-screen prompts to set up users.



Access Control Setup password screen



Access Control Setup screen

Notes

- Initial factory supervisor PIN = 5625
- The supervisor PIN cannot be deleted, and should be changed to prevent unauthorized user ID setup. The supervisor PIN does not allow access to the unit. At least one user ID must be set up to gain access to the unit.



Access Control keypad

Enter a valid PIN using the keypad.

5.2 Using i.Lock™ Interlock (*Pharmacy Pass-Thru only*)


The i.Lock™ Interlock feature is an electronically-activated magnetic lock which prohibits doors on opposite sides of the unit from being opened at the same time. This feature minimizes the risk of contaminating product or allowing airborne particulates to transfer from the uncontrolled side (ante-room or “dirty” room) to the clean room side of the unit.

Notes

- The interlock feature does not restrict the use of both the left and right doors on the same side of a double door unit.
- During a power failure, the optional door interlock feature will continue to prevent opposite doors from being opened simultaneously until battery power is depleted or until the back-up battery switch is switched OFF.
- During a power failure, switch the battery back-up switch OFF and use the mechanical door key to provide secure storage for refrigerator contents.
- Switching the back-up battery key switch OFF will disable the monitoring system during a power failure.


6 Horizon Series™ Operation

6.1 Initial Power-Up

1. Plug the power cord into a grounded outlet that meets the electrical requirements on the product specification label.
2. Switch the AC ON/OFF switch **ON**.
3. Install 9 V back-up battery provided (undercounter battery located in literature box; upright battery located on top of unit).
4. Switch the Alarm ON/OFF key switch to **ON**.
5. Press **Down Arrow** (Mute) if high temperature alarm sounds. 

Notes

- For models equipped with optional Access Control, switch the back-up battery key switch **ON**.
- During a power failure, the back-up battery continues to provide power to the optional Access Control lock (if equipped). If the back-up battery is not functioning, the optional Access Control lock will not secure the door.
- If an alarm condition other than High Temperature occurs, refer to the service manual for troubleshooting.

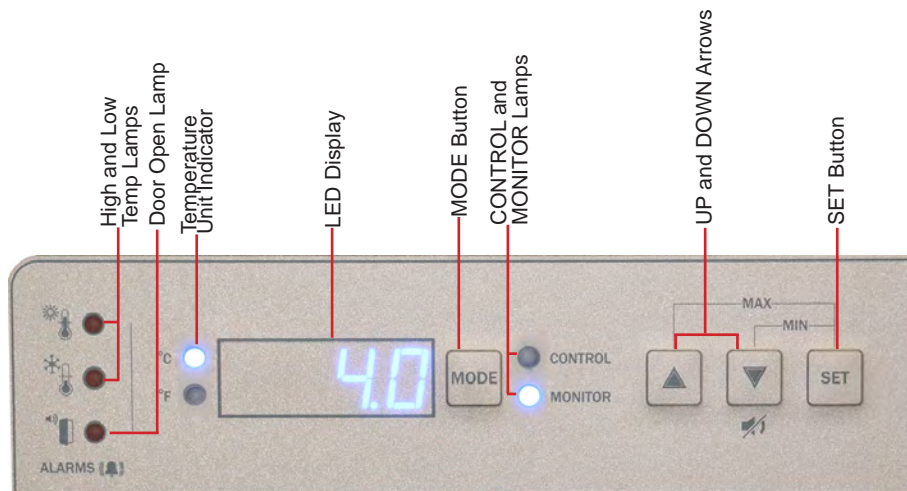


Horizon Series™ temperature monitor and control interface


6.2 Display Minimum and Maximum Monitor Temperature Recordings

The minimum and maximum recording feature allows the user to view a minimum temperature occurrence and a maximum temperature occurrence within a given period of time. The timer provides a time reference in which those temperatures occurred.

View minimum temperature recording

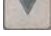1. Press and hold the **Down Arrow** button for 1 second and listen for a single beep. 
2. The display will alternate between **LO** and a valid temperature value five (5) times followed by a single beep to indicate exit back to the temperature display.



View maximum temperature recording

1. Press and hold the **Up Arrow** button for 1 second and listen for a single beep. 
2. The display will alternate between **HI** and a valid temperature value five (5) times followed by a single beep to indicate exit back to the temperature display.



View recorded temperature timer

Notes

- The timer denotes the period of time that has elapsed. It does not display the time at which a minimum or maximum temperature occurred.
- The maximum period of time the timer can record is 99:59 (99 hours and 59 minutes).

1. Press and hold either the **Up** or **Down Arrow** button for 1 second.  or .
2. While the display is flashing the **HI** or **LO** value, press and hold the **SET** button for 1 second.
3. The display will alternate five (5) times between **CLr** and a value representing the number of hours and minutes that have elapsed since the last recording (example: 12:47 would represent 12 hours and 47 minutes). A single beep will follow to indicate exit back to temperature display.

Clear minimum and maximum temperature recordings

1. Press and hold either the **Up** or **Down Arrow** button for 1 second.  or .
2. While the display is flashing the **HI** or **LO** value, press and hold the **SET** button for 1 second and listen for a single beep.
3. While the display is flashing the elapsed time since last reset, press and hold the **SET** button for 2 seconds. **CLr** will be displayed followed by a series of 3 beeps to indicate exit back to the temperature display.

Notes

The minimum and maximum temperature and timer will reset when:

- the unit is powered off and battery back-up is not engaged, or
- after 99 hours and 59 minutes have elapsed.

6.3 Change Temperature Setpoint

Note

Default setpoint is 4.0 °C for HBR and HLR models; 5.0 °C for HPR models.

1. Press and release **MODE** to change to Control mode. The CONTROL lamp will illuminate.
2. Press and hold **SET** to display the current setpoint temperature.
3. Hold **SET** and press the **Up** or **Down Arrow** as necessary to set the desired setpoint value.
4. Release all buttons; the setpoint is changed.
5. Press and release **MODE** to return to Monitor mode. The MONITOR lamp will illuminate.

6.4 Set Parameter Values

1. Press and hold the **Up** and **Down Arrows** simultaneously for 3 seconds to enter program mode.
2. The LED Display will show .C or .F to indicate Celsius or Fahrenheit.
3. Press and release **MODE** button to scroll through the parameters.
4. Once the desired parameter is selected, press and hold the **SET** button while pressing the **Up** or **Down Arrow** to select the desired value.
5. Release **SET** button. The new setting is saved.
6. Press and hold the **Up** and **Down Arrows** simultaneously for 3 seconds to exit program mode.

 **Note**

Contact Helmer Technical Service to set Rail Limit values.

Table 3. Parameter Values

Parameter	Visual Indicator	Range	Default
Celsius or Fahrenheit	None	°C, °F	°C
High Temperature	MONITOR Lamp & HIGH Lamp	-40.0 to 25.0 (°C) -40 to 77 (°F)	5.5 °C (<i>HBR and HLR models</i>)
			6.5 °C (<i>HPR models</i>)
Low Temperature	MONITOR Lamp & LOW Lamp	-40.0 to 25.0 (°C) -40 to 77 (°F)	1.5 °C (<i>HBR models</i>)
			2.0 °C (<i>HLR and HPR models</i>)
Monitor Offset	MONITOR Lamp	-10.0 to 10.0 (°C) -18 to 18 (°F)	Varies
Control Offset	CONTROL Lamp	-10.0 to 10.0 (°C) -18 to 18 (°F)	Varies
Upper Rail Limit	CONTROL Lamp & HIGH Lamp	0.1 to 10.0 (°C); 1 to 18 (°F)	0.7 °C
			1.5 °C (<i>Pass-Thru models</i>)
Lower Rail Limit	CONTROL Lamp & LOW Lamp	0.1 to 10.0 (°C); 1 to 18 (°F)	-0.7 °C
			-1.5 °C (<i>Pass-Thru models</i>)

6.5 Set Temperature Units

 **Note**

If temperature units are changed, the temperature setpoints, offsets and alarm settings must be recalibrated.

1. Press and hold the **Up** and **Down Arrows** simultaneously for 3 seconds to enter program mode.
2. The LED Display will show °C or °F to indicate Celsius or Fahrenheit.
3. Press and hold the **SET** button while pressing the **Up** or **Down Arrow** to select the desired temperature unit.
4. Release **SET** button. The new setting is saved.
5. Press and hold the **Up** and **Down Arrows** simultaneously for 3 seconds to exit program mode.

6.6 Set Alarm Setpoints (Parameters)

1. Press and hold the **Up** and **Down Arrows** simultaneously for 3 seconds to enter program mode.
2. The LED Display will show .C or .F to indicate Celsius or Fahrenheit.
3. Press **MODE** until HIGH TEMP or LOW TEMP and MONITOR lamps flash.
4. Hold **SET**, then press **Up** or **Down Arrow** to change the setpoint.
5. Release **SET** button. The new setting is saved.
6. Press and hold **Up** and **Down Arrows** simultaneously for 3 seconds to exit program mode.

6.7 Temperature Calibration Offsets

Temperature calibration offsets indicate an acceptable margin of error between the actual temperature value and the desired temperature value.

Monitor Offset

- ◆ Value is factory-set to match a calibrated reference thermometer.
- ◆ Refer to the service manual for instructions regarding changing the Monitor Offset.

Control Sensor Offset and Hysteresis

The control sensor affects the reading of the control probe temperature and therefore the actual temperature of the refrigerator. This should not be adjusted from the original setting unless directed by Helmer Technical Service.

The Upper and Lower Rail Limits help control the refrigeration based on the control probe temperature reading and the set point. These limit values should not be changed from the default setting unless directed by Helmer Technical Service.

NOTICE

Control Sensor Offset and Rail Limits are factory-preset and should not be changed. Contact Helmer Technical Service for instructions regarding changing these values.

6.8 Active Alarms

The controller displays temperature and alarm information.

Table 4. Horizon Series Active Alarms

Alarm	Visual Indicator	Description
High Temperature	HIGH TEMP lamp flashes	Chamber temperature reading is above high temperature alarm setpoint
Low Temperature	LOW TEMP lamp flashes	Chamber temperature reading is below low temperature alarm setpoint
Display/Control Board Communication Error	Er04	Display board fails to communicate with the control board
Control Board to Compressor Inverter Error	Er05	Communication loss from control board to compressor inverter
Power Failure	"PoFF" appears on display	Power to unit has been disrupted
Primary Monitor Probe Failure (RTD2)	Er02	Probe not functioning properly
Control Probe Failure (RTD1)	Er01	Probe not functioning properly
No Battery	bAtt	Battery voltage is low
Configuration Error	Er07	Indicates that an EEPROM reading was corrupted or dip switch settings on the control board have changed since last power-up
Door Open < 3 min.	DOOR ALARM lamp lights	Door is open (less than three minutes)
Door Open > 3 min.	DOOR ALARM lamp flashes	Door has been open 3 minutes or longer*

*Audible alarm will sound after door is open for 3 minutes.

6.9 Mute and Disable Audible Alarms

Note

Muting audible alarms does not disable alarm lamps or signals sent through the remote alarm interface.

- ◆ Press **Down Arrow** (Mute) to mute audible alarms.
- ◆ To disable all audible alarms, insert the key in the Alarm Disable switch and turn.

6.10 Light Operation (if installed)

The light switch is located on the monitoring and control panel and controls the LED light within the chamber.

7 Horizon Series™ Access Control (Optional)

Allows user-specific secure access to the refrigerator.

Notes

- During a power failure, the optional Access Control lock will remain locked until battery power is depleted or until the back-up battery key switch is switched OFF.
- During a power failure, switch the battery back-up switch OFF and use the mechanical door key to provide secure storage for refrigerator contents.
- Refer to Horizon Series Access Control manual for complete information.

7.1 Using Access Control

The Access Control keypad was programmed at the factory with a master code (0000). The master code is used to program the keypad and enter user codes.

Note

The master code should be changed to prevent unauthorized user code setup.

Enter unique user codes for up to 100 users. Each user code is stored with a specific record location number. Keep a log of the location numbers and user codes with users' names.

Add User Code

1. Enter the master code followed by the * (asterisk) key
2. Press **1** to initiate user code programming function
3. Enter the location number (00 - 99)
4. Enter the user code (4 - 8 digit number) followed by the # (pound) key
5. Press * (asterisk) to save changes and return to normal operation

Delete User Code

1. Enter the master code followed by the * (asterisk) key
2. Press **1** to initiate delete user code programming function
3. Enter the location number (00 - 99) followed by the # (pound) key
4. Press * (asterisk) to save changes and return to normal operation

Open Refrigerator with Access Control



1. Enter the user code
2. Press # (pound) key

8 Product Specifications

8.1 Operating Standards

These units are designed to operate under the following environmental conditions:

- ◆ Indoor use only
- ◆ Altitude (maximum): 2000 m (120,125, 245, and 256 models); 3000 m (105 and 113 models)
- ◆ Ambient temperature range: 15 °C to 32 °C (59°F to 90°F)
- ◆ Relative humidity (maximum for ambient temperature): 80% for temperatures up to 31 °C; decreasing linearly to 76% at 32 °C
- ◆ Temperature control range: 2 °C to 10 °C (35°F to 50°F)
- ◆ Overvoltage Category II
- ◆ Pollution Degree 2
- ◆ RF Emissions: Group 1 - Class A
- ◆ EMC Environment: Basic
- ◆ Sound level is less than or equal to 56 dB(A)

Table 5. Electrical Specifications (Laboratory, Blood Bank, and Pharmacy)

Model	Input Voltage & Frequency	Voltage Tolerance	Circuit Breakers	Current Draw	Power Source	Remote Alarm Capacity
105	115V 60 Hz	±10%	4A quantity 2	1.4A	Grounded outlet, meeting national electric code (NEC) in the U.S. and local electrical requirements in all locations.	115V or 230V: 1A at 33V (AC) RMS or 30V (DC)
	220-240V 50/60 Hz			0.85A		
113	115V 60 Hz		4A quantity 2	2.3A		
	220-240V 50/60 Hz			1.12A		
120	115V 60 Hz		7A quantity 2	2.8A		
	220-240V 50/60 Hz			1.55A		
125	115V 60 Hz		7A quantity 2	2.8A		
	220-240V 50/60 Hz			1.55A		
245	115V 60 Hz		7A quantity 2	4.3A		
	220-240V 50/60 Hz			2.5A		
256	115V 60 Hz		7A quantity 2	4.3A		
	220-240V 50/60 Hz			2.5A		
226	115V 60 Hz		9A quantity 2	4.6A		
	220-240V 50/60 Hz			3.1A		
458	115V 60 Hz	9A quantity 2	6.3A			
	220-240V 50/60 Hz		3.5A			

* Amperage values are subject to change. Refer to the product specification label on your unit for current values.

Notes

- The interface on the remote alarm monitoring system is intended for connection to the end user's central alarm system(s) that uses normally-open or normally-closed dry contacts.
- If an external power supply exceeding 33V (RMS) or 30V (DC) is connected to the remote alarm monitoring system's circuit, the remote alarm will not function properly or may be damaged.

8.2 Refrigerator Specifications

Table 6. Refrigerator Dimensions

Model	Voltage Code	Amps	Cu. Ft/ Liters	Cabinet	Door	Dimensions W x H x D in. (mm)	Net Wt. lbs (kg)
						Exterior*	
iLR105-GX	115V 60 Hz	1.4	5.3 (150)	Undercounter	Single hinged solid	24.3 x 31.9 x 28.1 (616 x 809 x 714)	174 (79)
	220-240V 50/60 Hz	0.85					
iBR105-GX	115V 60 Hz	1.4	5.3 (150)	Undercounter	Single hinged solid	24.3 x 31.9 x 28.1 (616 x 809 x 714)	186 (85)
	220-240V 50/60 Hz	0.85					
iPR105-GX	115V 60 Hz	1.4	5.3 (150)	Undercounter	Single hinged solid	24.3 x 31.9 x 28.1 (616 x 809 x 714)	182 (83)
	220-240V 50/60 Hz	0.85					
HLR105-GX	115V 60 Hz	1.4	5.3 (150)	Undercounter	Single hinged solid	24.3 x 31.9 x 27.8 (616 x 809 x 705)	174 (79)
	220-240V 50/60 Hz	0.85					
HBR105-GX	115V 60 Hz	1.4	5.3 (150)	Undercounter	Single hinged solid	24.3 x 31.9 x 27.8 (616 x 809 x 705)	184 (84)
	220-240V 50/60 Hz	0.85					
HPR105-GX	115V 60 Hz	1.4	5.3 (150)	Undercounter	Single hinged solid	24.3 x 31.9 x 27.8 (616 x 809 x 705)	182 (83)
	220-240V 50/60 Hz	0.85					
iLR113-GX	115V 60 Hz	2.3	13 (377)	Slimline	Single hinged glass	24.6 x 70.5 x 31.0 (625 x 1790 x 787)	306 (139)
	220-240V 50/60 Hz	1.12					
iBR113-GX	115V 60 Hz	2.3	13 (377)	Slimline	Single hinged glass	24.6 x 70.5 x 31.0 (625 x 1790 x 787)	342 (156)
	220-240V 50/60 Hz	1.12					
iPR113-GX	115V 60 Hz	2.3	13 (377)	Slimline	Single hinged glass	24.6 x 70.5 x 31.0 (625 x 1790 x 787)	338 (154)
	220-240V 50/60 Hz	1.12					
HLR113-GX	115V 60 Hz	2.3	13 (377)	Slimline	Single hinged glass	24.6 x 70.5 x 31.0 (625 x 1790 x 787)	306 (139)
	220-240V 50/60 Hz	1.12					
HBR113-GX	115V 60 Hz	2.3	13 (377)	Slimline	Single hinged glass	24.6 x 70.5 x 31.0 (625 x 1790 x 787)	337 (153)
	220-240V 50/60 Hz	1.12					
HPR113-GX	115V 60 Hz	2.3	13 (377)	Slimline	Single hinged glass	24.6 x 70.5 x 31.0 (625 x 1790 x 787)	338 (154)
	220-240V 50/60 Hz	1.12					
iLR120-GX	115V 60 Hz	2.8	20 (572)	Upright	Single hinged glass	29.5 x 79.6 x 31.9 (748 x 2021 x 810)	445 (202)
	220-240V 50/60 Hz	1.55					
iBR120-GX	115V 60 Hz	2.8	20 (572)	Upright	Single hinged glass	29.5 x 79.6 x 31.9 (748 x 2021 x 810)	507 (230)
	220-240V 50/60 Hz	1.55					
iPR120-GX	115V 60 Hz	2.8	20 (572)	Upright	Single hinged glass	29.5 x 79.6 x 31.9 (748 x 2021 x 810)	438 (199)
	220-240V 50/60 Hz	1.55					
HLR120-GX	115V 60 Hz	2.8	20 (572)	Upright	Single hinged glass	29.5 x 78.3 x 31.9 (748 x 1989 x 810)	442 (201)
	220-240V 50/60 Hz	1.55					
HBR120-GX	115V 60 Hz	2.8	20 (572)	Upright	Single hinged glass	29.5 x 78.3 x 31.9 (748 x 1989 x 810)	504 (229)
	220-240V 50/60 Hz	1.55					
HPR120-GX	115V 60 Hz	2.8	20 (572)	Upright	Single hinged glass	29.5 x 78.3 x 31.9 (748 x 1989 x 810)	487 (221)
	220-240V 50/60 Hz	1.55					
iLR125-GX	115V 60 Hz	2.8	25 (714)	Upright	Single hinged glass	29.5 x 79.6 x 37.9 (748 x 2021 x 962)	456 (207)
	220-240V 50/60 Hz	1.55					
iBR125-GX	115V 60 Hz	2.8	25 (714)	Upright	Single hinged glass	29.5 x 79.6 x 37.9 (748 x 2021 x 962)	535 (243)
	220-240V 50/60 Hz	1.55					
iPR125-GX	115V 60 Hz	2.8	25 (714)	Upright	Single hinged glass	29.5 x 79.6 x 37.9 (748 x 2021 x 962)	517 (235)
	220-240V 50/60 Hz	1.55					
HLR125-GX	115V 60 Hz	2.8	25 (714)	Upright	Single hinged glass	29.5 x 78.3 x 37.9 (748 x 1989 x 962)	453 (206)
	220-240V 50/60 Hz	1.55					
HBR125-GX	115V 60 Hz	2.8	25 (714)	Upright	Single hinged glass	29.5 x 78.3 x 37.9 (748 x 1989 x 962)	532 (242)
	220-240V 50/60 Hz	1.55					

Model	Voltage Code	Amps	Cu. Ft/ Liters	Cabinet	Door	Dimensions W x H x D in. (mm)	Net Wt. lbs (kg)
						Exterior*	
HPR125-GX	115V 60 Hz	2.8	25 (714)	Upright	Single hinged glass	29.5 x 78.3 x 37.9 (748 x 1989 x 962)	514 (234)
	220-240V 50/60 Hz	1.55					
iLR245-GX	115V 60 Hz	4.3	45 (1271)	Upright	Double hinged glass	59.0 x 79.6 x 31.9 (1499 x 2021 x 810)	667 (303)
	220-240V 50/60 Hz	2.5					
iBR245-GX	115V 60 Hz	4.3	45 (1271)	Upright	Double hinged glass	59.0 x 79.6 x 31.9 (1499 x 2021 x 810)	809 (367)
	220-240V 50/60 Hz	2.5					
iPR245-GX	115V 60 Hz	4.3	45 (1271)	Upright	Double hinged glass	59.0 x 79.6 x 31.9 (1499 x 2021 x 810)	775 (352)
	220-240V 50/60 Hz	2.5					
HLR245-GX	115V 60 Hz	4.3	45 (1271)	Upright	Double hinged glass	59.0 x 78.3 x 31.9 (1499 x 1989 x 810)	667 (303)
	220-240V 50/60 Hz	2.5					
HBR245-GX	115V 60 Hz	4.3	45 (1271)	Upright	Double hinged glass	59.0 x 78.3 x 31.9 (1499 x 1989 x 810)	808 (367)
	220-240V 50/60 Hz	2.5					
HPR245-GX	115V 60 Hz	4.3	45 (1271)	Upright	Double hinged glass	59.0 x 78.3 x 31.9 (1499 x 1989 x 810)	774 (352)
	220-240V 50/60 Hz	2.5					
iLR256-GX	115V 60 Hz	4.3	56 (1586)	Upright	Double hinged glass	59.0 x 79.6 x 37.9 (1499 x 2021 x 962)	703 (319)
	220-240V 50/60 Hz	2.5					
iBR256-GX	115V 60 Hz	4.3	56 (1586)	Upright	Double hinged glass	59.0 x 79.6 x 37.9 (1499 x 2021 x 962)	863 (392)
	220-240V 50/60 Hz	2.5					
iPR256-GX	115V 60 Hz	4.3	56 (1586)	Upright	Double hinged glass	59.0 x 79.6 x 37.9 (1499 x 2021 x 962)	827 (376)
	220-240V 50/60 Hz	2.5					
HLR256-GX	115V 60 Hz	4.3	56 (1586)	Upright	Double hinged glass	59.0 x 78.3 x 37.9 (1499 x 1989 x 962)	693 (315)
	220-240V 50/60 Hz	2.5					
HBR256-GX	115V 60 Hz	4.3	56 (1586)	Upright	Double hinged glass	59.0 x 78.3 x 37.9 (1499 x 1989 x 962)	853 (387)
	220-240V 50/60 Hz	2.5					
HPR256-GX	115V 60 Hz	4.3	56 (1586)	Upright	Double hinged glass	59.0 x 78.3 x 37.9 (1499 x 1989 x 962)	817 (371)
	220-240V 50/60 Hz	2.5					
iBR226-GX	115V 60 Hz	4.6	26.4 (747)	Upright	Single hinged glass x 2	29.9 x 79.6 x 40.9 759 x 2021 x 1038	548 (249)
	220-240V 50/60 Hz	3.1					
iPR226-GX	115V 60 Hz	4.6	26.4 (747)	Upright	Single hinged glass x 2	29.9 x 79.6 x 40.9 759 x 2021 x 1038	526 (239)
	220-240V 50/60 Hz	3.1					
HBR226-GX	115V 60 Hz	4.6	26.4 (747)	Upright	Single hinged glass x 2	29.1 x 79.3 x 40.0 740 x 2013 x 1016	541 (246)
	220-240V 50/60 Hz	3.1					
HPR226-GX	115V 60 Hz	4.6	26.4 (747)	Upright	Single hinged glass x 2	29.1 x 79.3 x 40.0 740 x 2013 x 1016	523 (238)
	220-240V 50/60 Hz	3.1					
iBR458-GX	115V 60 Hz	6.3	58.5 (1656)	Upright	Double hinged glass x 4	59 x 79.6 x 40.9 1499 x 2021 x 1038	913 (414)
	220-240V 50/60 Hz	3.5					
iPR458-GX	115V 60 Hz	6.3	58.5 (1656)	Upright	Double hinged glass x 4	59 x 79.6 x 40.9 1499 x 2021 x 1038	868 (394)
	220-240V 50/60 Hz	3.5					
HBR458-GX	115V 60 Hz	6.3	58.5 (1656)	Upright	Double hinged glass x 4	59 x 79.3 x 40.0 1499 x 2013 x 1016	899 (408)
	220-240V 50/60 Hz	3.5					
HPR458-GX	115V 60 Hz	6.3	58.5 (1656)	Upright	Double hinged glass x 4	59 x 79.3 x 40.0 1499 x 2013 x 1016	862 (391)
	220-240V 50/60 Hz	3.5					

* 105 models - Exterior dimensions include handle and leveling feet at lowest level.

** 113, 120, 125, 245, 256 models - Exterior dimensions include casters, door handle and electrical box.

*** 226, 458 models - Exterior dimensions include casters, door handles (all models), and i.Series bezel (i.Series models only).

Table 7. Storage Component Specifications

Storage Component	Net weight lbs (kg) <i>105 models</i>	Net weight lbs (kg) <i>113 models</i>	Net weight lbs (kg) <i>120/245 models</i>	Net weight lbs (kg) <i>125/256 models</i>	Net weight lbs (kg) <i>226/458 models</i>
Shelf	6 (3)	6 (3)	7 (3.2)	8 (3.6)	8.7 (4)
Ventilated Drawer	10 (5)	8.2 (3.7)	11 (5)	13.5 (6.1)	-
Liquid-tight Stainless Drawer (<i>i.Series</i>)	12 (6)	12.1 (5.5)	14.2 (6.5)	17.2 (7.8)	16.1 (7.3)
Liquid-tight Stainless Drawer (<i>Horizon Series</i>)	11 (5)	11.8 (5.4)	13.7 (6.2)	16.6 (7.5)	15 (6.8)
Stainless Drawer with Locking Lid	15 (7)	13.6 (6.2)	15.5 (7)	19 (8.6)	-

 **Notes**

- Amperage values listed represent the highest current draw presented among available factory configurations for each model Units without heated glass doors will have lower current draw.
- Maximum height added with leveling feet or casters installed is 2" (51 mm).
- Maximum load per shelf or drawer - 100 lbs (46kg).
- Net weight may vary depending on storage configuration. Weight listed in the table reflects standard configuration for each model.

9 Compliance

9.1 Safety Compliance



This device complies with the requirements of directive 93/42/EEC concerning Medical Devices, as amended by 2007/47/EEC 2017/745 for Medical Devices Article 120.



This product is certified to applicable UL and CSA standards by a NRTL.

This product is IECEE CB Scheme certified and complies with national differences for safety certification beyond IEC 61010-1-12 3rd edition.

9.2 Environmental Compliance



This device complies with the 2011/65/EU Directive for the Restriction of Hazardous Substances (RoHS) Recast including the amendment to Annex II described in Commission Delegated Directive (EU) 2015/863.



This device falls under the scope of Directive 2012/19/EU Waste Electrical and Electronic Equipment (WEEE) .

When disposing of this product in countries affected by this directive:

- ◆ Do not dispose of this product as unsorted municipal waste.
- ◆ Collect this product separately.
- ◆ Use the collection and return systems available locally.

For more information on the return, recovery, or recycling of this product, contact your local distributor.

9.3 EMC Compliance

Helmer Scientific Refrigerators meet the applicable requirements of IEC61326 and EN55011 and are intended for use in the electromagnetic environment specified in the Operating Standards section of this manual. The customer or the user of these devices should assure they are used in such environment.



This device complies with FCC Radiated and Conducted Emissions Approval to CFR47, Part 15; Class A levels

Helmer Scientific

14400 Bergen Boulevard, Noblesville, IN 46060 USA

Life Science Solutions

A business unit of Trane Technologies