

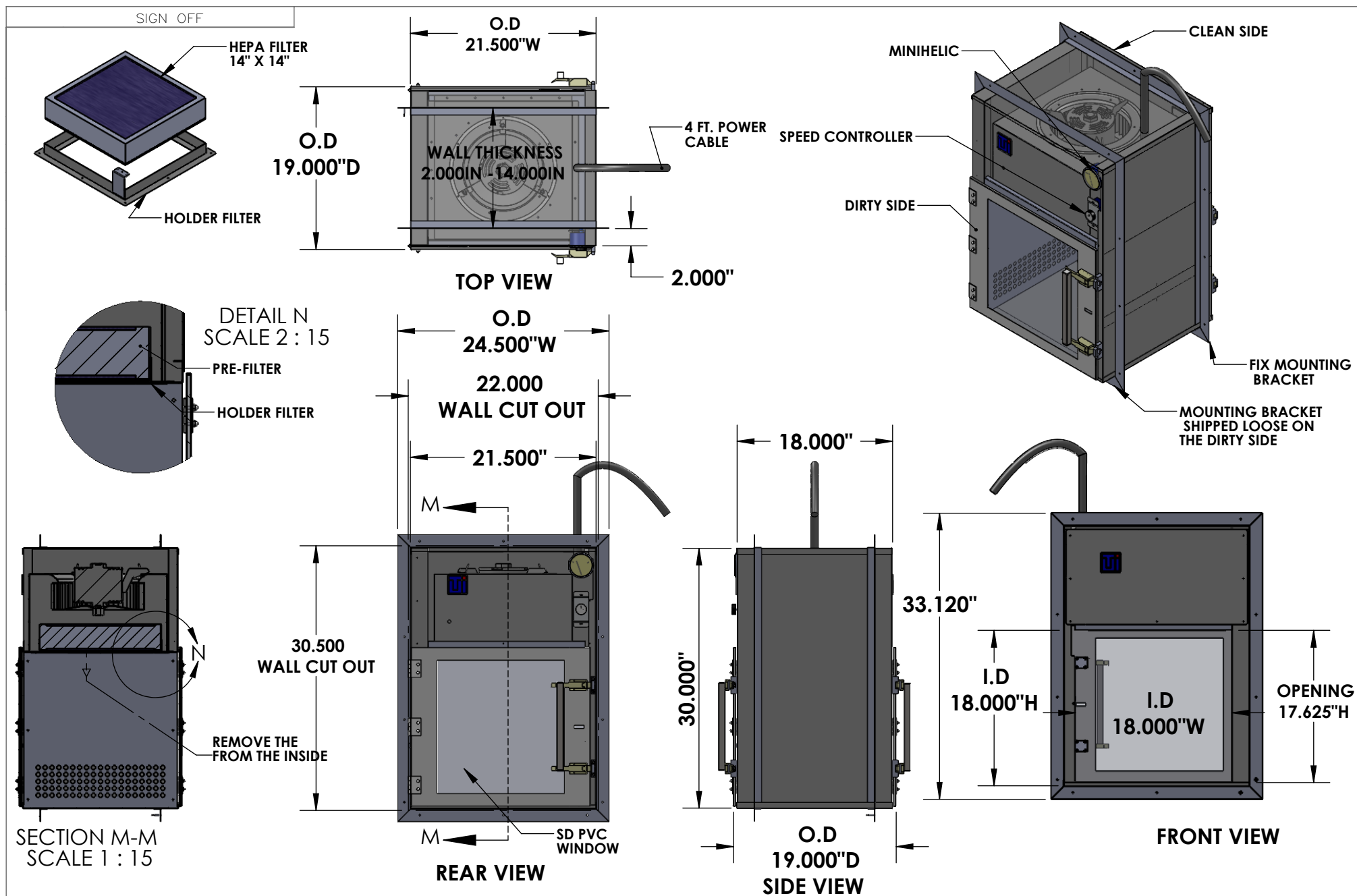


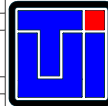


REVISIONS				MATERIAL SPEC. 304SS	FINISH SPEC.	UNLESS OTHERWISE SPECIFIED, DIMENSIONS ARE IN INCHES		 Terra Universal PROPRIETARY & CONFIDENTIAL ALL RIGHTS RESERVED BY TERRA UNIVERSAL DISCLOSURE OR COPYING PROHIBITED	TITLE		
RV#	DESCRIPTION	DATE	BY			TOLERANCES ON:			A 2634-80-A	PASS-THROUGH; RECIRCULATING HEPA FILTRATION, 18"W X 18"D X 18"H WALL MOUNT, 304 SS	
NC	INITIAL RELEASE	10-9-19	JR			.XX DECIMALS ±.03	X" ± 1"				
						.XXX DECIMALS ±.010	X/X ± 1/64				
						SIGNATURES					
						DATE					
						ISSUED BY: JUAN REAL					
						10-9-19					
				APPROVED BY:							
				REF #: 255754-N1-1							
								SCALE		SHEET 1 OF 1	


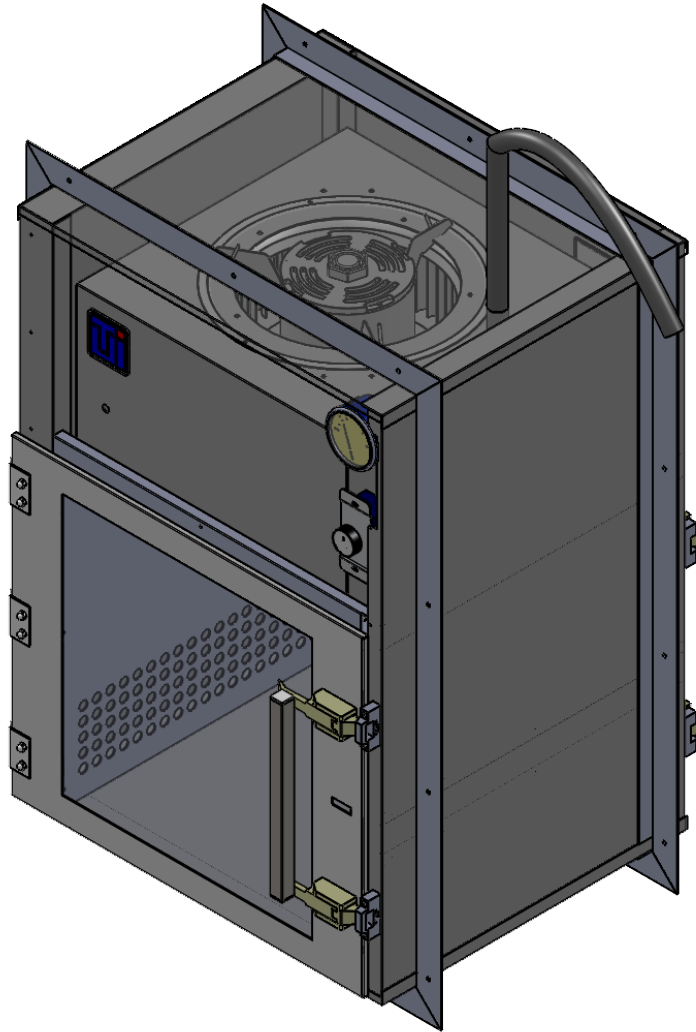


Recirculating Pass-Through Specifications

Electrical Specifications	120 VAC 1 Amp 1 PH
Chamber Dimension	18"W X 18"D X 18"H
Wall Cutout	22.000"W x 30.500"H
Overall Dimension	21.5"W X 20.5"D X 30"H
Rough Opening Dimension	18"W X 18"H
Window Dimension	14.75"W X 15.375"H
Window Material	Static Dissipative PVC
Blower Flow Rate	105 CFM
HEPA Filter Dimension	14"W X 14"D X 3.25"H
HEPA Filter Specifications	99.999 Efficiency @ 0.3 Microns
Material	304 Stainless Steel 18 Gauge
Fire Rating	None
Door Gasket Material	One-Piece "e" profile Santoprene non-outgassing rubber gasket with fused corners
HEPA Filter Replacement Location	Inside the Pass-Through


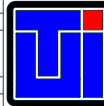
REVISIONS				MATERIAL SPEC.	FINISH SPEC.	UNLESS OTHERWISE SPECIFIED, DIMENSIONS ARE IN INCHES				 Terra Universal PROPRIETARY & CONFIDENTIAL ALL RIGHTS RESERVED BY TERRA UNIVERSAL DISCLOSURE OR COPYING PROHIBITED	TITLE	
RV#	DESCRIPTION	DATE	BY			TOLERANCES ON:					SCALE	SHEET 2 OF 4
NC	INITIAL RELEASE	8/30/19	JP			.XX DECIMALS ±.03 X' ± 1' .XXX DECIMALS ±.010 X/X' ± 1/64						
						SIGNATURES		DATE				
						ISSUED BY: JUAN REAL		8/30/19				
						APPROVED BY:						
						REF #:		253647				
											PASS-THROUGH; RECIRCULATING HEPA FILTRATION, 18"W X 18"D X 18"H WALL MOUNT, 304 SS	
											2634-80-3016566	


STEP 1:
DRILL A HOLE IN THE
WALL FOR THE POWER CABLE.
THEN SLIDE THE PASS-THROUGH
INTO THE WALL AND SCREW IT ON
THE ATTACHED MOUNTING BRACKET
TO HOLD IT IN PLACE

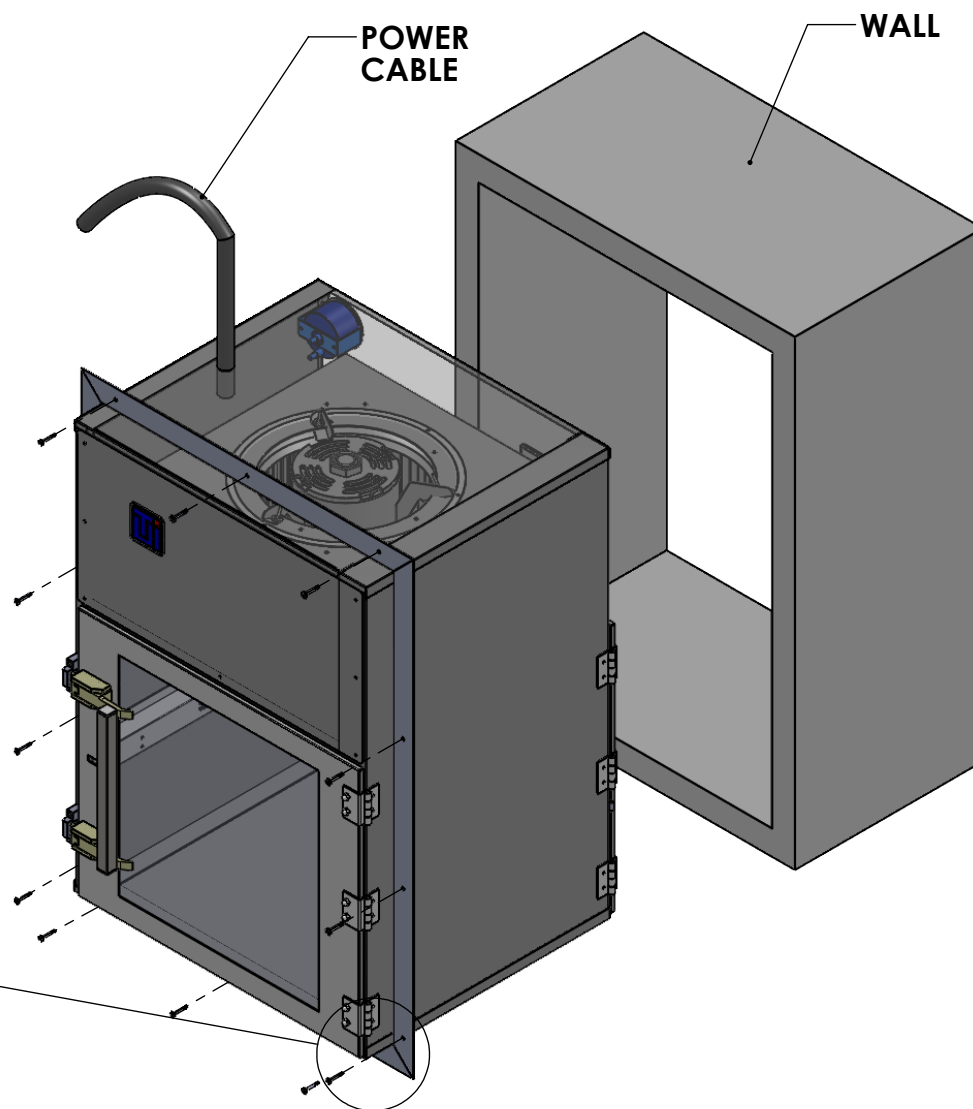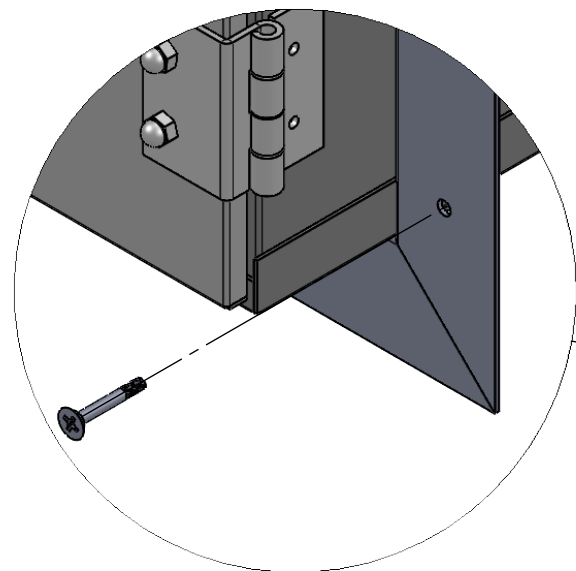

REVISIONS				MATERIAL SPEC. 304SS	FINISH SPEC.	UNLESS OTHERWISE SPECIFIED, DIMENSIONS ARE IN INCHES				 Terra Universal PROPRIETARY & CONFIDENTIAL ALL RIGHTS RESERVED BY TERRA UNIVERSAL DISCLOSURE OR COPYING PROHIBITED	TITLE					
RV#	DESCRIPTION	DATE	BY			TOLERANCES ON:					PASS-THROUGH; RECIRCULATING HEPA FILTRATION, 18"W X 18"D X 18"H WALL MOUNT, 304 SS					
NC	INITIAL RELEASE	8-30-19	JR			.XX DECIMALS ±.03 X" ± 1" .XXX DECIMALS ±.010 X/X ± 1/64										
						SIGNATURES					DATE					
						ISSUED BY: JUAN REAL					8-30-19					
						APPROVED BY:										
						REF #:					253647 EC-1					
											A	2634-80-3016566				
											SCALE				SHEET 3 OF 4	

STEP 2:
SLIDE THE MOUNTING BRACKET
INTO THE WALL ON THE DIRTY SIDE
AND SCREW IT ON TO HOLD IT IN PLACE

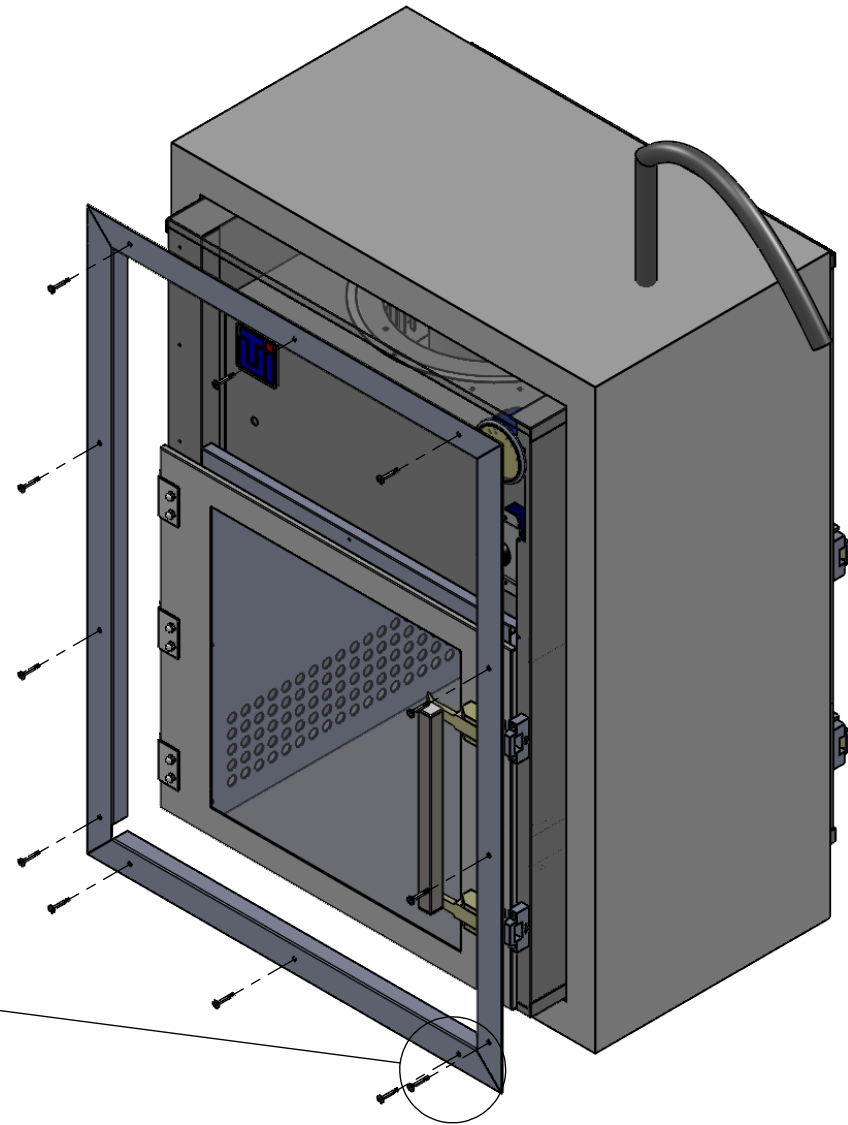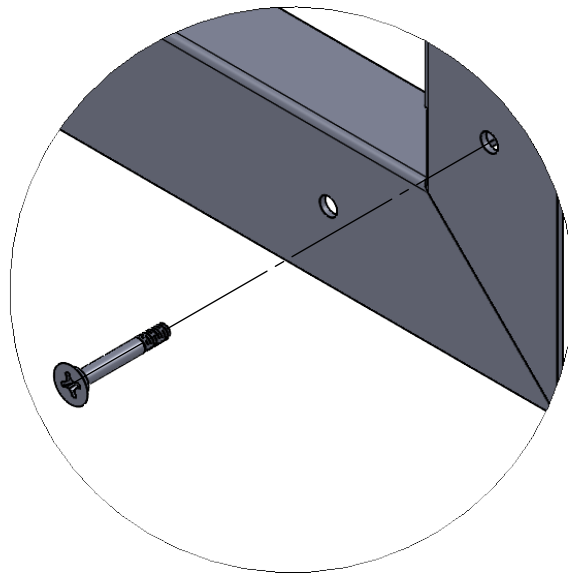

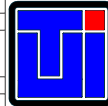REVISIONS				MATERIAL SPEC. 304SS	FINISH SPEC.	UNLESS OTHERWISE SPECIFIED, DIMENSIONS ARE IN INCHES				 Terra Universal PROPRIETARY & CONFIDENTIAL ALL RIGHTS RESERVED BY TERRA UNIVERSAL DISCLOSURE OR COPYING PROHIBITED	TITLE		
RV#	DESCRIPTION	DATE	BY			TOLERANCES ON:					PASS-THROUGH; RECIRCULATING HEPA FILTRATION, 18"W X 18"D X 18"H WALL MOUNT, 304 SS		
NC	INITIAL RELEASE	8/30/19	JP										
						.XX DECIMALS ±.03 X' ± 1'							
						.XXX DECIMALS ±.010 X/X ± 1/64							
						SIGNATURES		DATE			2634-80-3016566		
						ISSUED BY: JUAN REAL		8/30/19					
						APPROVED BY:					A		
				REF #:	253647				SCALE		SHEET 4 OF 4		