



Helmer Scientific, Noblesville, Indiana
www.helmerinc.com

Technical Data Sheet

High-Performance, Medical-grade Refrigerator


Model: HLR125-GXC, Horizon Series™ Laboratory / Upright

Application, Rating, & Electrical Data	
Application	High performance storage of medical and scientific products
Storage Volume	25.2 cu ft 714 L
Temperature Range	+2°C to +10°C
Set Point	+4°C
Electrical Power	115V 60Hz 208-240V 50/60Hz
Maximum Current	2.8A 1.55A
Facility Supply Rating	Branch circuit with earth ground, meeting plug/cord ratings as well as local electrical safety codes and requirements 110-120V AC: 15 A (minimum) 208-240V AC: 10 A (minimum)
Power Plug/Power Cord Length	NEMA 5-15 hospital-grade 115V 60Hz, 10 ft (3.0m) Inquire about additional plug/cord options
Certification/Agency Listing	QPS (Certified to UL and CSA Standards) IEC/UL61010-2-011:2016 IEC 61010-1:2010, AMD1:2016
Energy Star	Yes, ENERGY STAR® Certified High Performance Refrigerator
Indoor/Outdoor Use	Indoor use only, non-residential
Application Environment	Non-corrosive, non-flammable, non-explosive
Ambient Operating Temperature	+15°C to +32°C +59°F to +90°F
Refrigeration	
Refrigeration System	VCC compressor-based, forced-air hydrocarbon
Compressor	Hermetic, variable speed (VCC). Rated speed range: 1300-4000 rpm
Condenser Type	Air-cooled fin and tube, sealed
Expansion Device	Cap tube
Evaporator Type	Direct expansion forced-air fin and tube
Defrost Method	Automatic
Refrigerant	R600a, Isobutane, EPA SNAP compliant
Controller / Configuration Settings	
User Interface	LED digital display, °C or °F
Power Switch	On/Off
Controller Type	Microprocessor with button input and LED display
Security	Lockable door
Control Sensor	NTC, sealed
Communication Ports	Remote alarm contacts
Power Failure Alarm	Yes
High/Low Alarms	Yes, fully adjustable
Door Ajar Alarm	Yes
Min/Max Temperature	Yes, display and reset
Temperature Monitoring Ballast	Yes, 4 oz bottle with screw cap, solid ballast optional
Battery Back-up	Yes, 9 hr monitoring system Rechargeable 12V, lead acid


Performance Characteristics	
Uniformity °C	+/-1.0
Recovery to 8°C after 3 min door opening min	3
Stability °C	0.05
Energy Consumption kWh/day	3.01
Noise emission dB	49
Heat Rejection / Emission BTU/hr	915
Pull-down time to 4°C min	43

All performance data from 25 cu ft refrigerator, 24.8°C ambient, 4.0°C


Dimensions and Construction	
Interior (w x h x d)	24.8 x 58.1 x 30 in 628 x 1476 x 762 mm
Exterior (w x h x d)	29.0 x 79.0 x 34.1 in 737 x 2006 x 866 mm
Overall Exterior (w x h x d)	29.0 x 79.0 x 38.5 in 737 x 2006 x 977 mm
Insulation	Sustainable, US EPA and SNAP approved foam
Exterior/Interior Finish	Bacteria-resistant powder-coating
Door	1, dual-pane glass, right-hinge
Access Port	Top Access Port- for external monitoring probe(s) (0.75 in 19 mm diameter) Left Sidewall Port - with interior and exterior plugs (1.75 in 45 mm diameter)
Interior Storage (w x d)	3 ventilated, powder coated shelves, adjustable 22.5 x 28 in 572 x 711 mm, 100 lb 46 kg max capacity/shelf
Casters	Yes, swivel-locking
Net Weight	400 lb 182 kg
Shipping Weight	462 lb 210 kg
Clearance Requirements	Minimum of 8 in (203 mm) above and 3 in (76 mm) behind unit for proper ventilation, clearance, and feature access
Options Accessories	Temperature Chart recorder (stand-alone), Chart paper, Leveling feet, Solid ballast, Remote alarms, Floor and Wall bracket kit, Seismic anchoring Kit, Storage trays, Left-hinge door, Solid door, Extended warranty
Warranty	Rel.i™ Warranty- 5 years compressor, 2 years parts, 1 year labor Outside US and Canada- Contact your distributor for warranty information

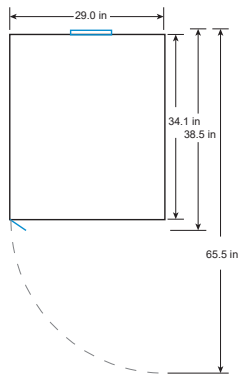


TOP VIEW

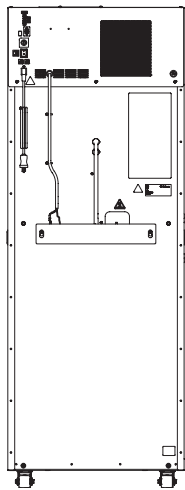

BACK VIEW

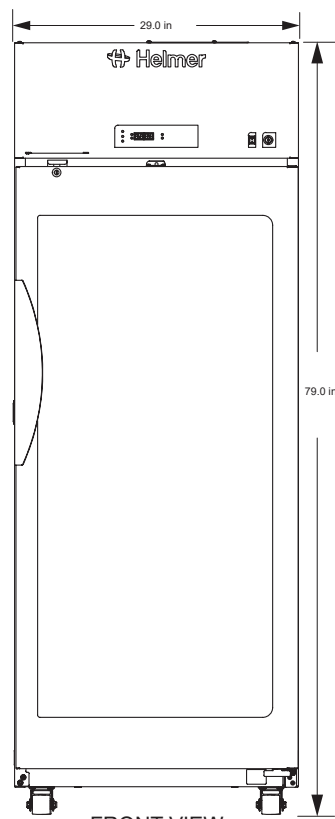

FRONT VIEW

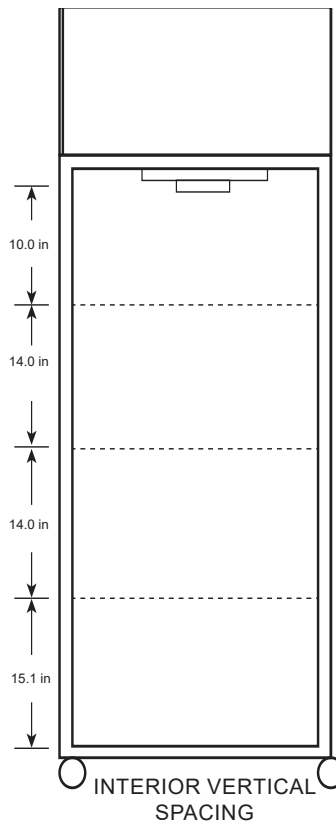

INTERIOR VERTICAL SPACING